HEINNER

CÂNTAR CORPORAL DE DIAGNOSTIC

Model: HDS-150BKSL

- Cântar corporal de diagnostic
- Capacitate: 150 kg

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Importator: **Network One Distribution**

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ELECTRONIC PERSONAL SCALE

Model: HDS-150BKS

- Electronic personal scale
- Capacity: 150kg

Thank you for purchasing this product!

I. INTRODUCTION

Before using the appliance, please read carefully this instruction manual and keep it for future reference.

This instruction manual is designed to provide you with all required instructions related to the installation, use and maintenance of the appliance.

In order to operate the unit correctly and safely, please read this instruction manual carefully before installation and usage.

II. CONTENT OF THE PACKAGE

- **Electronic personal scale**
- **Instruction manual**
- **Warranty card**

III. PRODUCT INTRODUCTION

ABOUT THE PRODUCT

The body mass index (BMI) is calculated from a person's weight and height and is a widely accepted measurement of weight. A BMI of 20~25 is considered to be ideal. Lower than 20 is too slim and higher than 25 is considered over weight. (BMI is only for reference).

ADVISE FOR USE & CARE

1. The glass platform will be slippery when it is wet. Keep it dry! Keep the scale in cool place with dry air.
2. Keep standing still during the course of weighing. DO NOT strike, shock or drop the scale.
3. Treat the scale with care as it is a precise instrument. Always keep your scale in horizontal position.
4. Clean the scale with a damp cloth and avoid water getting inside the scale. DO NOT use chemical/abrasive cleaning agents.
5. If failed to turn on the scale, check whether battery is installed, or battery power is low. Insert or replace new battery.
6. If there is error on display or cannot turn off the scale for a long time, please take off the battery for about 3 seconds, then install it again to remove soft-fault. If you cannot resolve the problem, please contact the seller for consultation.
7. Do not use it for commercial purpose.

IV. USING INSTRUCTIONS

LCD DISPLAY


KEY ILLUSTRATION

5. (Guest): press to enter guest setting mode when the scale is off/shows zero.
6. (Set): when the scale is off/shows zero, press to set user parameters including user code, gender, height; when the scale shows zero, long push the button to switch the unit.
7. (▲): advance the value by one interval each press. Continuous pressing will speed up the advance of the value. Press to enter the user test mode when weighing.
8. (▼): decrease the value by one interval each press. Continuous pressing will speed up the decrease of the value. Press to enter the user test mode when weighing.

WEIGHING

Auto step on weighing: Add on weighing object (>3kg) to turn on the scale, it then starts to weigh and shows the total weight. The LCD flashes three times after stabilizes, and search saved users' data automatically. If the difference between current weight and saved user data is smaller than 5kg, it will calculate the BMI automatically with saved user data. If more than one user suits this condition, the user code will flash then you can choose your own user code. If no action, it will go back to the weighing mode and auto off after 10 seconds.

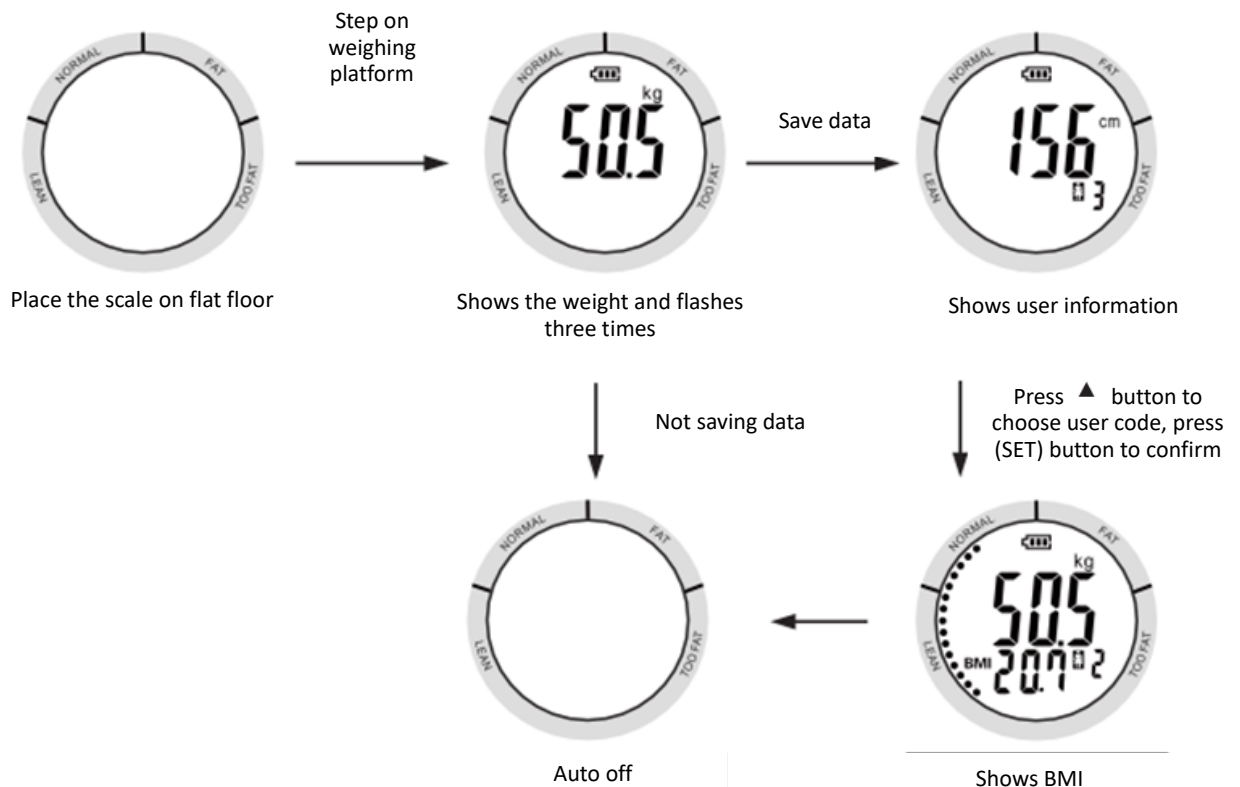After digits stabilize, more additional loadings can be conducted. Ensure each add-on loading has to exceed 2kg to get a new reading.

Example: if you want to get a weighing value of your baby, you can weigh yourself first and then hold your baby together, without stepping off the scale. Then you will get the total weight reading. Subtracting the previous reading of yourself from the new reading, you will get the difference as your baby's weighing value.


BMI MODE

Parameter can be set in the following range:

User code	gender	height
1~8	male/female	100~250 cm 3'03.0"~8'02.5"


2. BMI user mode (save data automatically)

To set user parameters: press SET to set user parameter, setting in order of user code (1-8), gender, height. Press SET again to confirm the input. When finish setting, press the scale gently, the parameters will be saved and the scale goes into weighing mode. Step on the scale, when digit stabilizes, it shows the weight and BMI and the circle bar will show the BMI level.

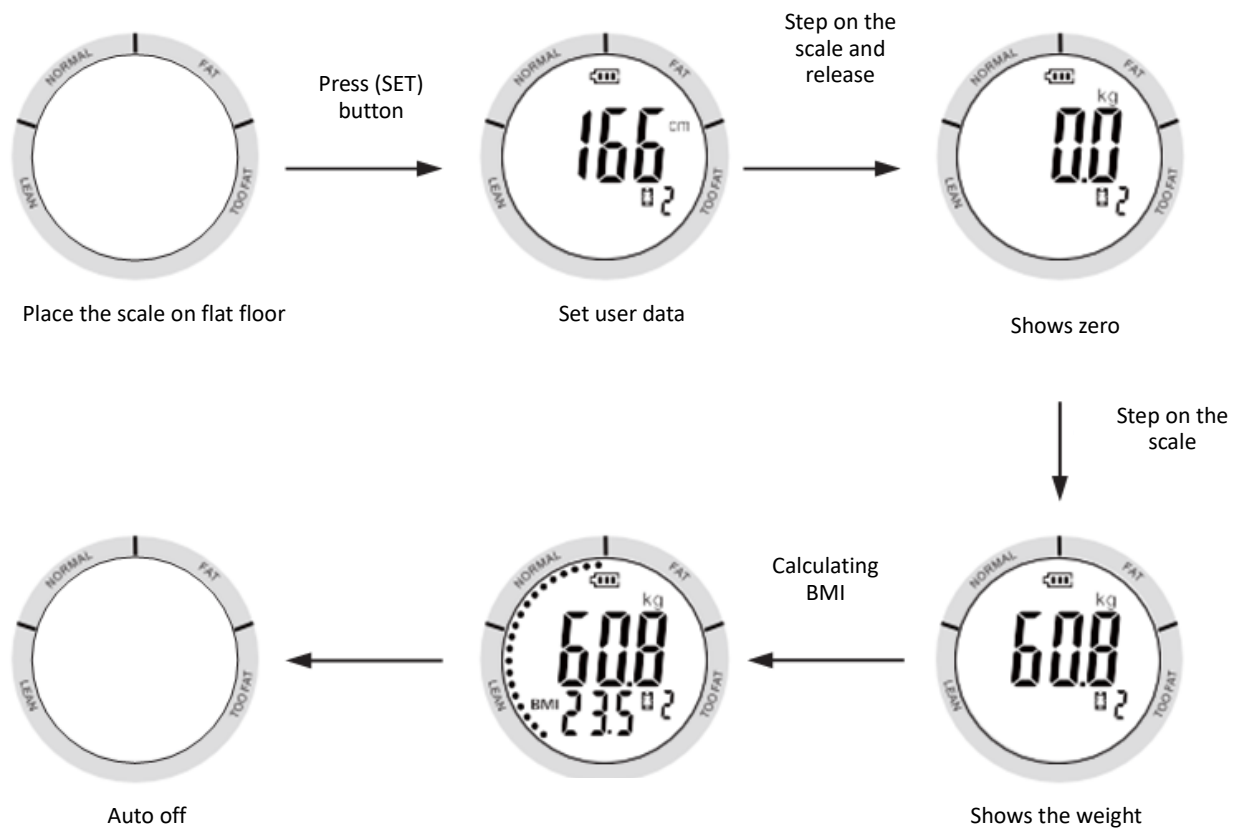
2. BMI guest mode (data will not be saved)

To set user parameters: Press GUEST button to turn on the scale and enter GUEST mode, the LCD shows "GUS". Setting gender and height, press set again to confirm. When finish setting, press the scale gently to start weighing. Step on the scale, when digit stabilizes, it shows the weight and BMI, and the circle bar will show the BMI level.

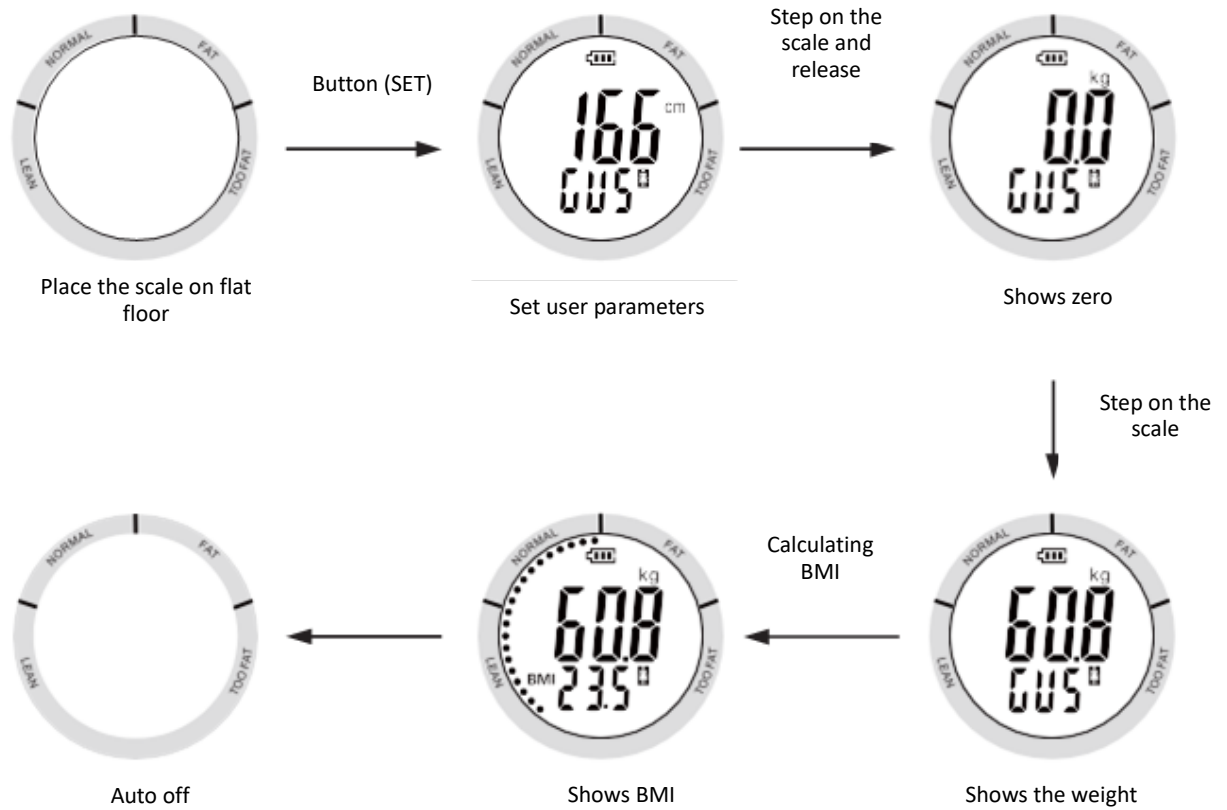3. BMI user test mode

- i: Under weighing mode, press ▲ button or ▼ button when shows zero, the scale chooses user code automatically and enter weight and BMI mode.
- ii: Step on the scale after selection, the scale shows weighing results when digit stabilizes.

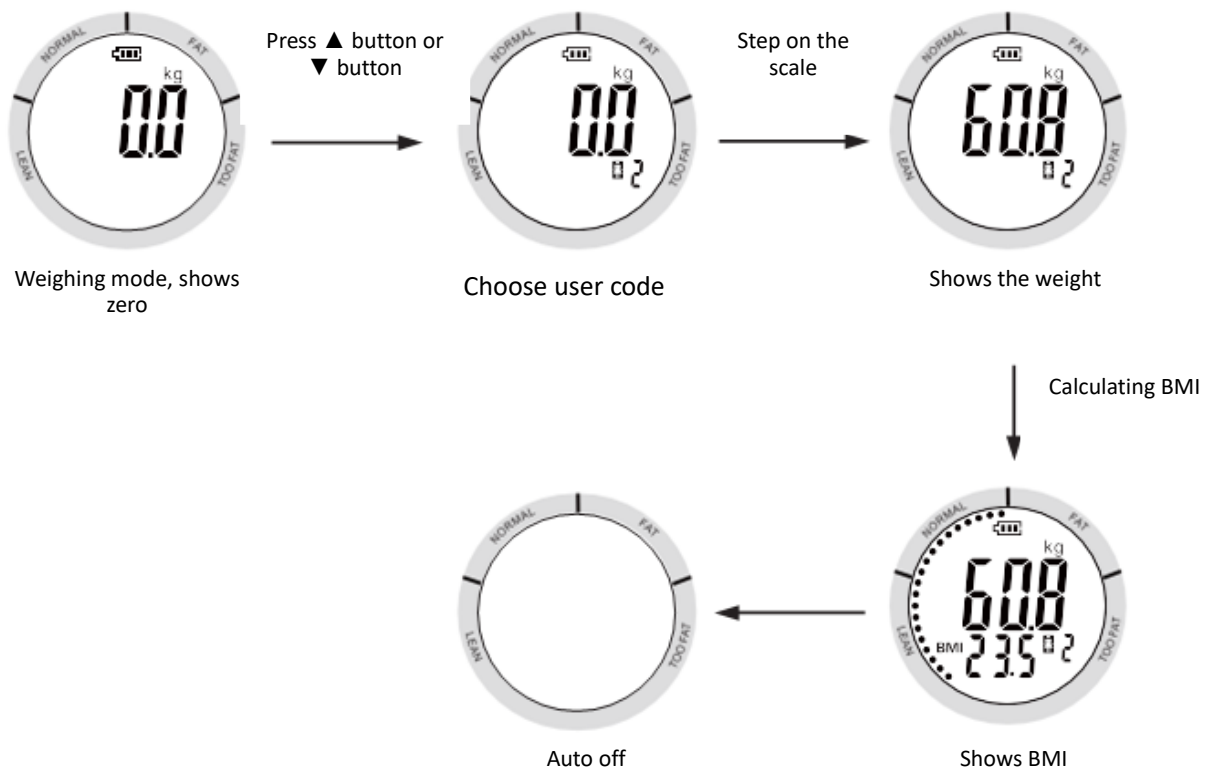WARNING INDICATORS

OVERLOAD INDICATION

Scale overload (max. weight capacity is indicated on the scale). Please step off to avoid damage.

LOW BATTERY INDICATION

Battery power is low. Please replace the battery.

RE-TEST


Display "L" when you step off the scale, means this weight reading with error. Please re-test.

BATTERY

The battery type of the scale is one of the following:

- c. 1 x 3V CR2032 lithium batteries. Remove isolation sheet before use;
- d. 2 x 1.5V AAA batteries. Remove package foil before use.


REPLACING THE BATTERIES

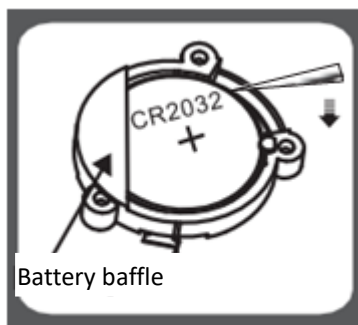PIC 1

Condition 1: open battery case cover on the back of scale, if battery case is the type shown in picture 1.

- C. Take out the used battery with the aid of a sharp object, as shown in picture 1.
- D. Install new battery by putting one side of the battery beneath the battery baffle and then pressing down the other side. Please take note to the polarity.

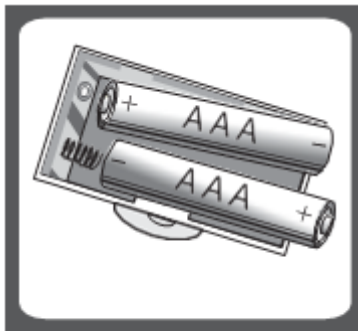PIC 2

Condition 2: Open battery case cover on the back of scale, if battery case is the type shown in picture 2.

- C. Take out the used battery (with the aid of a sharp object, if needed)
- D. Install new battery by putting one side of the battery and then pressing down the other side. Please take note to the polarity.

V. TECHNICAL FICHE

Capacity	150 kg
Battery	2 x 1.5V AAA
Color	Black/silver

**Environment friendly disposal**

You can help protect the environment!

Please remember to respect the local regulations: hand in the non-working electrical equipment's to

HEINNER is a registered trademark of **Network One Distribution SRL**. Other brands and product names are trademarks or registered trademarks of their respective holders.

No part of the specifications may be reproduced in any form or by any means or used to make any derivative such as translation, transformation, or adaptation without permission from NETWORK ONE DISTRIBUTION.

Copyright © 2013 Network One Distribution. All rights reserved.

www.heinner.com, <http://www.nod.ro>
<http://www.heinner.com>, <http://www.nod.ro>

This product is in conformity with norms and standards of European Community

Importer: **Network One Distribution**

Str. Marcel Iancu, nr. 3-5, Bucuresti, Romania

Tel: +40 21 211 18 56, www.heinner.ro, www.nod.ro

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Copyright © 2013 Network One Distribution. Toate drepturile rezervate.